

ECHOSAMPLER

INTRACAVITAL PATHOLOGY AND ENDOMETRIAL HISTOLOGY

The EchoSampler is a flexible sampling catheter used for dual diagnosis: intra cavital pathology and endometrial histology.

FEATURES


- dual diagnosis: cavital pathology and endometrial histology
- for ultrasound guided biopsy
- for saline infusion sonohysterography
- implemented in case of abnormal uterine bleeding
- flexible sampling catheter
- depth markings for ease in locating depth of insertion in the uterus
- tissue collection and saline infusion through side hole
- smooth rounded distal tip and small size for minimal patient discomfort and easy passage through cervical canal
- handle with luer lock for a more secure fitting to syringes
- echogenic line to facilitate visualisation when using ultrasound



ORIGINAL PICTURE

GYNO INFO ECHOSAMPLER

Length	180 mm (± 3 mm)
Tube OD	2.93 mm (± 0.04 mm)
Sterile	ETO sterilized
Packaging	75/shipping box
Order number	#5164 ☎

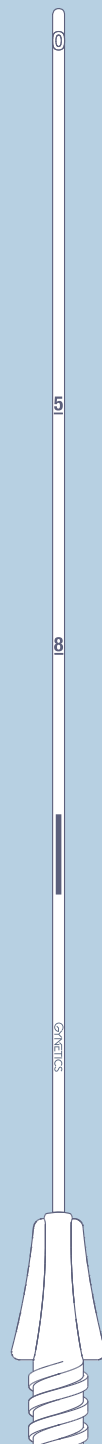

GMP Belgium Echosampler 02/2017