


Oosafe®

OOSAFE® AIR CLEANER KN

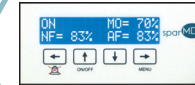

sparMED

GOosafe in IVF Laboratories

# OOSAFE® AIR CLEANER KN

OOSAFE® AIR CLEANER KN IS ADVANCED EQUIPMENT, DESIGNED FOR IVF LABORATORY USE, IT PURIFIES AIR AGAINST VOLATILE ORGANIC COMPOUNDS (VOCs), CHEMICALLY ACTIVE COMPOUNDS (CACs), PARTICLES, MICROORGANISMS, TOXINS, SOLVENTS AND ODOR.



Oosafe® Air Cleaner KN  
(including initial filter set)  
Product code:  
**OOACEQ-02**



Manufactured according  
to ISO 9001 Quality  
Management System

## OOSAFE® AIR CLEANER KN

- **Better air quality, better embryo development:** Oosafe® Air Cleaner KN purifies the air in laboratories and improves embryo development. Oocytes and embryos do not have an immune system to fight against toxic compounds.
- **More activated carbon, more VOC removal capacity:** Canister Carbon Filter Set consists of 5 pieces of activated carbon canisters made of pure activated carbon which gives unmatched VOCs removal capacity.
- **No chemicals inside:** High quality purified air is guaranteed by pure materials only. Chemicals inside the filters may be dangerous to human health and embryo development.
- **Pure activated carbon, the only way to trap VOCs:** VOCs cannot be trapped by HEPA filters because their size is smaller than the HEPA filter's pore size. VOCs can only be trapped inside activated carbon.
- **45° air blowing:** Faster air purification with 45° air blowing system.
- **Large HEPA Filter:** Large HEPA filter for longer replacement period and better air quality, purifies air against particles and microorganisms.
- **Alarm function:** Specific sensors are used to define the exact time when filters should be replaced.
- **Low noise level:** Higher user comfort is achieved with very low noise level (lower than 50dB (A) at full speed and 42dB (A) at half speed).
- **Approved by scientists:** Oosafe® Air Cleaner KN has been used for many years in well-known IVF clinics all over the world as the main air purification device.
- **Initial filter set:** Oosafe® Air Cleaner KN comes with Main Prefilter, Canister Carbon Filter Set and HEPA filter.
- **Optional:** 4 stage filtering with Dual Extra Prefilter is possible.

# OOSAFE® AIR CLEANER KN

- **Adjustable and Accurate Air Speed Level:**  
Air speed can be adjusted from the LCD Screen.
- **No Money Wasted on Replacement of Filters:**  
Use of highly sensitive air speed sensor plus timer alerts means filters are replaced when necessary - no more replacement of filters before they are filled.
- **Easy to operate:**  
Main service can be handled on the LCD screen through service menu. Replacement of the filters can be done by the users.

## ORDER INFORMATION

**OOACEQ-02** Oosafe® Air Cleaner KN including initial filter set  
**OOACHF-02** Oosafe® HEPA Filter for Oosafe® Air Cleaner KN  
**OOACCF-01** Oosafe® Canister Carbon Filter Set for Oosafe® Air Cleaner KN  
**OOACPF-01** Oosafe® Main Prefilter for Oosafe® Air Cleaner KN  
**OOACPF-02** Oosafe® Extra Prefilter for Oosafe® Air Cleaner KN (Optional)



Oosafe® Air Cleaner KN  
Product code:  
**OOACEQ-02**

TECHNICAL DATA	Dimensions	660x533x1640 (w x d x h, in mm)
	Air exchanges per hour (AEH)	9-28 times
	Weight (including filters)	77 kg
	HEPA Filter efficiency (MPPS)	≥ 99.995%
	Total amount of activated carbon	9.4 kg


## INITIAL FILTERS



Oosafe® HEPA Filter for  
Oosafe® Air Cleaner KN  
Product code:  
**OOACHF-02**




Oosafe® Canister Carbon Filter Set for  
Oosafe® Air Cleaner KN  
Product code:  
**OOACCF-01**

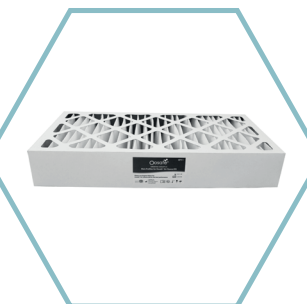

Oosafe® Main Prefilter for  
Oosafe® Air Cleaner KN  
Product code:  
**OOACPF-01**

## OPTIONAL

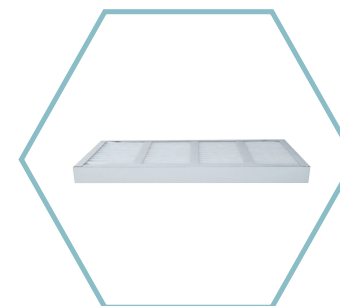

Oosafe® Extra Prefilter for  
Oosafe® Air Cleaner KN  
Product code:  
**OOACPF-02**



SparMED ApS - Ryttermarken 2, 3520 Farum, Denmark - Phone (+45) 39 40 25 03 - [www.sparmed.dk](http://www.sparmed.dk)



 [facebook.com/sparMED](https://facebook.com/sparMED)

 [twitter.com/sparMED](https://twitter.com/sparMED)

 [linkedin.com/company/sparmed-aps](https://linkedin.com/company/sparmed-aps)